

# Drive with 501RL2 pumphead



300 series cased pump

## Features and benefits

- Manual or remote 4-20 mA or 0-10 V analogue auto-control or RS232 control commands
- Up to 133:1 enhanced digital manual speed control
- Keypad lock prevents tampering or accidental changes
- MemoDose facility for accurate single-shot dosing
- Digital TTL inputs for run/stop and direction
- Instantly reversible
- Two-year warranty



## Technical specifications

|                                  | Drive with 501RL2 pumphead                                      |
|----------------------------------|---|
| Pumphead number of channels      | 1   |
| Speed control ratio              | 147:1   |
| Max. operating speed             | 200 rpm   |
| Min. operating speed             | 1.5 rpm   |
| Operating temperature range      | 5 to 40 °C  |
| Operating temperature range      | 40 to 104 °F  |
| Weight                           | 5.6 kg  |
| Weight                           | 12.3 lbs  |
| Drive control options            | U   |
| Drive standards                  | CE, cETLus  |
| Drive ingress rating             | IP31  |
| Drive humidity                   | 80% up to 31°C, 88°F, decreasing linearly to 50% at 40°C, 104°F |
| Drive noise                      | <70dB(A) at 1m  |
| Max. altitude                    | 2000  |
| Drive power supply               | 1ph, 100-120 V, 220-240 V 50/60 Hz, 100 VA                      |
| Compatible tubing bore size      | 1.6, 3.2, 4.8, 6.4, 8 mm  |
| Compatible tubing wall thickness | 2.4 mm  |


## Materials of construction

|                         | Drive with 501RL2 pumphead   |
|-------------------------|--|
| Drive casework          | ABS (Acrylonitrile butadiene styrene), Powder coated aluminum LM24 |
| Guide roller assembly   | Delrin   |
| Pumphead guard          | Polycarbonate  |
| Pumphead rotor assembly | Mazak, chromium plated   |
| Spindles                | Stainless steel  |
| Springs                 | Stainless steel  |
| Tube clamp assembly     | Acetal copolymer (Kematal)   |

Information listed covers the complete range.

For detailed specifications of individual models/components refer to user manual or contact WMFTS representative.

## Drive with 501RL2 pumphead dimensions

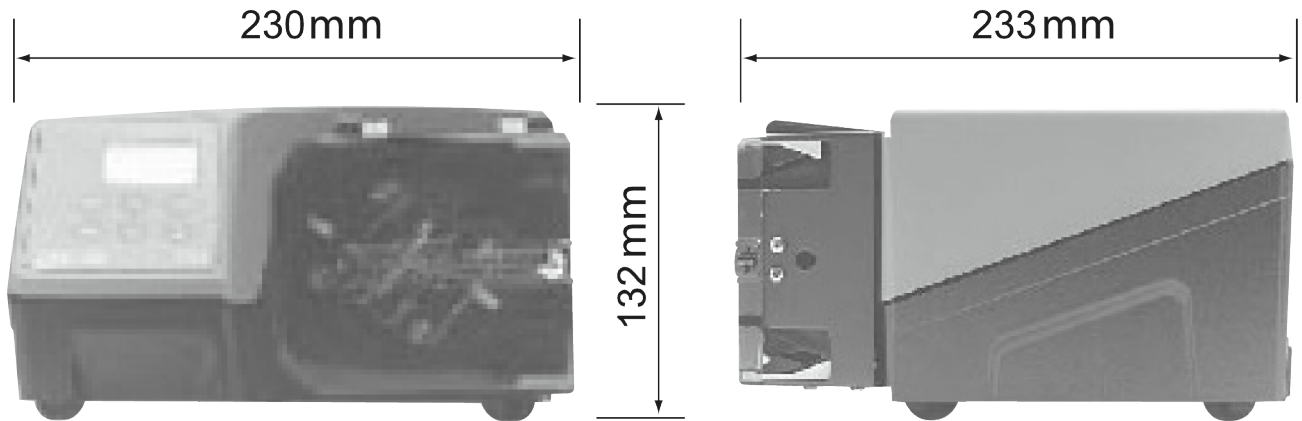

## Control options

| Functionality  | 323E | 323S | 323U | 323Du |
|--|------|------|------|-------|
| <b>Manual control</b>                                    |      |      |      |       |
| Keypad speed control                                     | ✓    | ✓    | ✓    | ✓     |
| 15-400rpm 27:1 speed control                             | ✓    |      |      |       |
| 3-400rpm 133:1 speed control                             |      | ✓    | ✓    | ✓     |
| 1.5-220rpm 147:1 speed control                           |      | ✓    | ✓    | ✓     |
| MemoDose for accurate single shot dosing                 |      | ✓    | ✓    | ✓     |
| <b>Remote control</b>                                    |      |      |      |       |
| Run/Stop, direction change via contact closure or 5V TTL |      |      | ✓    | ✓     |
| Remote operation of MemoDose (hand/foot switch)          |      |      | ✓    | ✓     |
| <b>Analogue speed control</b>                            |      |      |      |       |
| Inputs; 0-10V or 4-20mA                                  |      |      | ✓    | ✓     |
| <b>RS232 digital communication</b>                       |      |      |      |       |
| Digital serial control                                   |      |      |      | ✓     |
| <b>Security</b>  |      |      |      |       |
| Keypad lock to protect set-up                            |      | ✓    | ✓    | ✓     |

## Product codes

| Pumpheads product codes                |              |
|--|--------------|
| Description                            | Partcode     |
| 501RL2   pumphead for 2.4 mm wall tube | 053.0001.L20 |

\* Replace last 0 with A, E or U for American, European or UK mains supply (American 110 V 1ph 60 Hz)

If required for use with STA-PURE and Chem-Sure tubing, please contact Watson-Marlow Bredel technical support for drive and pumphead selection

| Tubing product codes |              |                     |                     |              |
|----------------------|--------------|---------------------|---------------------|--------------|
| Bore/wall (mm)       | Bioprene     | Marprene            | Platinum Silicone   |              |
| 1.6 / 2.4            | 903.0016.024 | 902.0016.024        | 913.A016.024        |              |
| 3.2 / 2.4            | 903.0032.024 | 902.0032.024        | 913.A032.024        |              |
| 4.8 / 2.4            | 903.0048.024 | 902.0048.024        | 913.A048.024        |              |
| 6.4 / 2.4            | 903.0064.024 | 902.0064.024        | 913.A064.024        |              |
| 8.0 / 2.4            | 903.0080.024 | 902.0080.024        | 913.A080.024        |              |
| Bore/wall (mm)       | STA-PURE®    | STA-PURE®Series PCS | STA-PURE®Series PFL | Chem-Sure®   |
| 1.6 / 2.4            | 960.0016.024 | 961.0016.024        | 966.0016.024        | 965.0016.024 |
| 3.2 / 2.4            | 960.0032.024 | 961.0032.024        | 966.0032.024        | 965.0032.024 |
| 4.8 / 2.4            | 960.0048.024 | 961.0048.024        | 966.0048.024        | 965.0048.024 |
| 6.4 / 2.4            | 960.0064.024 | 961.0064.024        | 966.0064.024        | 965.0064.024 |
| 8.0 / 2.4            | 960.0080.024 |                     |                     | 965.0080.024 |

---

Disclaimer: All flow rates shown were obtained pumping water at 20 °C (68 °F) with zero suction and delivery heads. The information contained in this document is believed to be correct but Watson-Marlow Limited accepts no liability for any errors it contains and reserves the right to alter specifications without notice. It is the users responsibility to ensure product suitability for use within their application. Watson-Marlow, LoadSure, Pumpsil, PureWeld XL, Bioprene, Marprene are registered trademarks of Watson-Marlow Limited. Tri-Clamp is a registered trademark of Alfa Laval Corporate AB. GORE and STA-PURE are registered trademarks of W.L. Gore and Associates. Please state the product code when ordering pumps and tubing.

[wmfts.com/global](https://wmfts.com/global)



07 August 2023