

Cellefill oRAB

Cellefill™

Features and benefits

- Fast track your cell and gene and biological therapies and futureproof for GMP compliance
- Maximise production efficiency with recipe driven, remote set-up and no format parts for the entire vial range
- Highly accurate, precision filling from micro volumes to 50 ml with better than $\pm 1\%$ achievable at most volumes
- Fully integrated environmental monitoring system in strategically risk-assessed locations
- Reduced risk of batch loss with flexible in-process controls including 100% weight check with zero waste start-up
- Integrated aseptic vial filling and containment solution reduces risk and cost
- Expert guidance and support throughout your project, from initial enquiry through to installation, to after sales support and validation package expertise, we are your project partner every step of the way.




Technical specifications

	Cellefill oRAB
Filling volume	Limited only by RTU vial size, filling volume range is from micro-volumes to 50 ml
Accuracy	achievable at most volumes $\pm 1.0\%$
Max Capacity	1500 units/hour
Batch Size	10 - 10,000 vials per batch
Operating temperature	5 °C to 40 °C (41 °F to 104 °F)
Power supply	230/400 V 3 ph, 50 Hz
Weight	2900 kg (6393 lbs)
Load per square meter	605 Kg
Barrier Technology	Open Restricted Access Barrier
Filling Line Type	Aseptic process filling line with manual disinfection
Application	Sterile product protection of sensitive biologics, antibodies, cell protein or nucleic acid-based products
Cleanroom Class	Grade B / ISO 7 cleanroom
Vial Size	ISO 2R - 50R
Weigh Check	20-100 %
Stoppers	ISO standard or Lyo, 13 mm and 20 mm
Infeed and outfeed trays	Vial RTU compatible with SG EZ-Fill® and SCHOTT adaptiQ® trays
Port type	Castus® alpha RTP 190 iCON, Sartorius Biosafe® 110 Monolever
Validation	Turnkey Validation FAT, SAT IQOQ

Throughput capacity up to 1500 vials/hour (based on batch size of 1000 vials/hour and dependent on customer configuration). NTT process for de-bagging inner bag (oRAB model) or both inner/out bags (iASP/iBIO versions) Environmental monitoring: Integrated viable and non-viable, active and passive Glove testers: optional accessory for all models Automated hygrometer control system: Optional on iASP and iBIO models

Cellefill oRAB dimensions

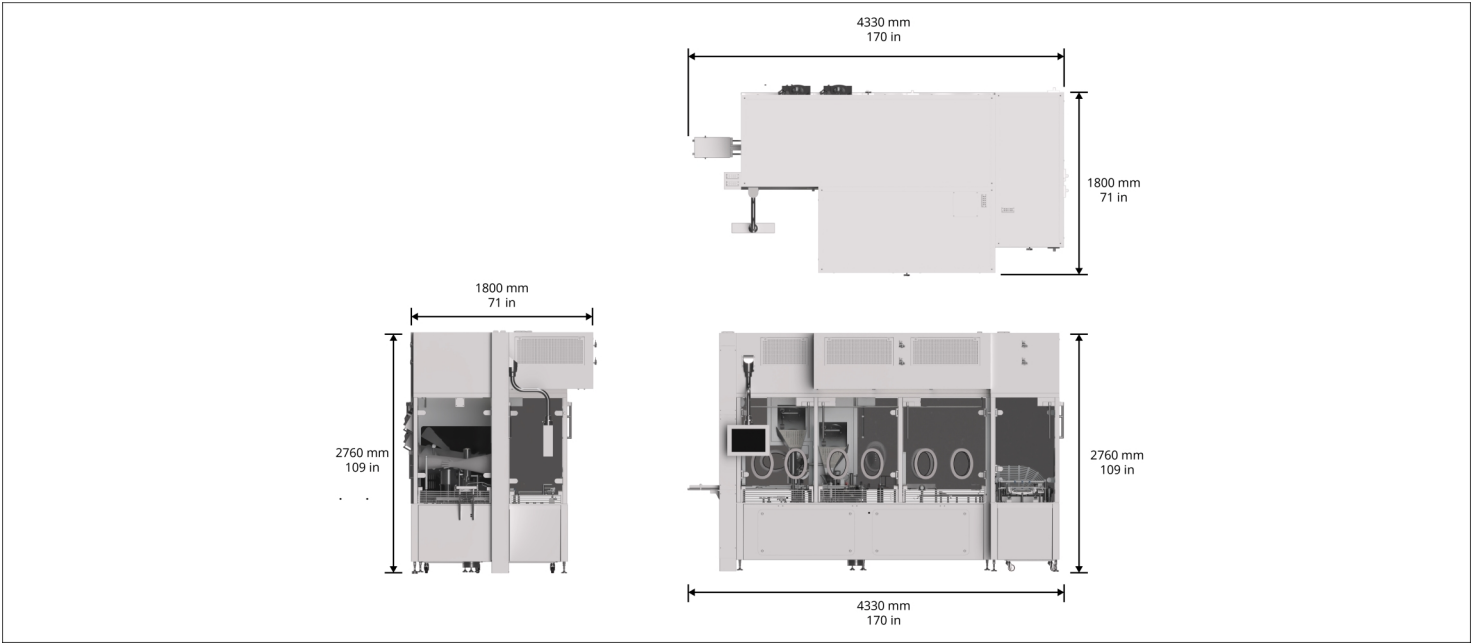

Description	oRAB
Minimum required operational space	2000 mm
Footprint LOA	4330 mm
Footprint WOA	1800 mm
Height Overall	2760 mm

Disclaimer: The information in this document is believed to be correct, but Watson-Marlow Flexicon A/S accepts no liability for any errors it contains, and reserves the right to alter specifications without notice. **WARNING:** These products are not designed for use in, and should not be used for, patient-connected applications. Cellefill is a registered trademark of Watson-Marlow Limited. SG EZ-Fill is a registered trademark of Stevanato Group. adaptiQ is a registered trademark of SCHOTT Pharma.



wmfts.com/global
17 December 2025