

Point-of-use valves

Point-of-use valves

Features and benefits

- Designed to promote drainage, eliminate dead-legs and process steam stagnation
- Provide a port to allow steam or material to be purged rather than trapped
- U-shaped valve design draining or sampling
- Integrated elbows for connection to piping loops
- Minimal L/D ratio enables easy cleaning
- Flush-mounted to the tank, no dead legs or dead spaces
- Easy to clean and self-draining by design - can be CIP/SIP to eliminate caking
- Available valve sizes: 0.5", 0.75", 1.0", 1.5", and 2.0"



Performance

Point-of-use valve flow rates		
Outlet Size (inches)	Cv at 1 psi (0.07 bar)	
	GPM	LPM
0.50	2.7	10.2
0.75	10.6	40
1.00	15.8	59.8
1.50	48	180
2.00	72	272


Technical specifications

	Point-of-use valves
Mounting type	Flush
Flow rate	10.2 - 272 l/m
Flow rate	2.7 - 72 USGPM
Sizes	0.5", 0.75", 1.0", 1.5", 2.0"
Outlet connections	Buttweld, Sanitary flange
Max. operating pressure	250 psi (17 bar)
Operating temperature	-60 °F to 350 °F (-51 °C to 177 °C)
Compatible actuators	Compact pneumatic, Fail closed, Fail open, Manual
Compatible diaphragms	EPDM, PTFE, Silicone, Silicone plus, Viton
Surface finish options	Electropolished, Max 10 micro-inch Ra (0.25 µm Ra), Max 15 micro-inch Ra (0.375 µm Ra), Max 20 micro-inch Ra (0.5 µm Ra)
Certification and compliance	ISO 9001
Standards	ASME BPE, CE-PED
Weight	3.33 - 37.43 lbs (1.5 - 17 kg)

Materials of construction

	Point-of-use valves
Material	Hastelloy C22, Hastelloy C276, Stainless steel 316L, Stainless steel AL-6XN

Point-of-use valves dimensions

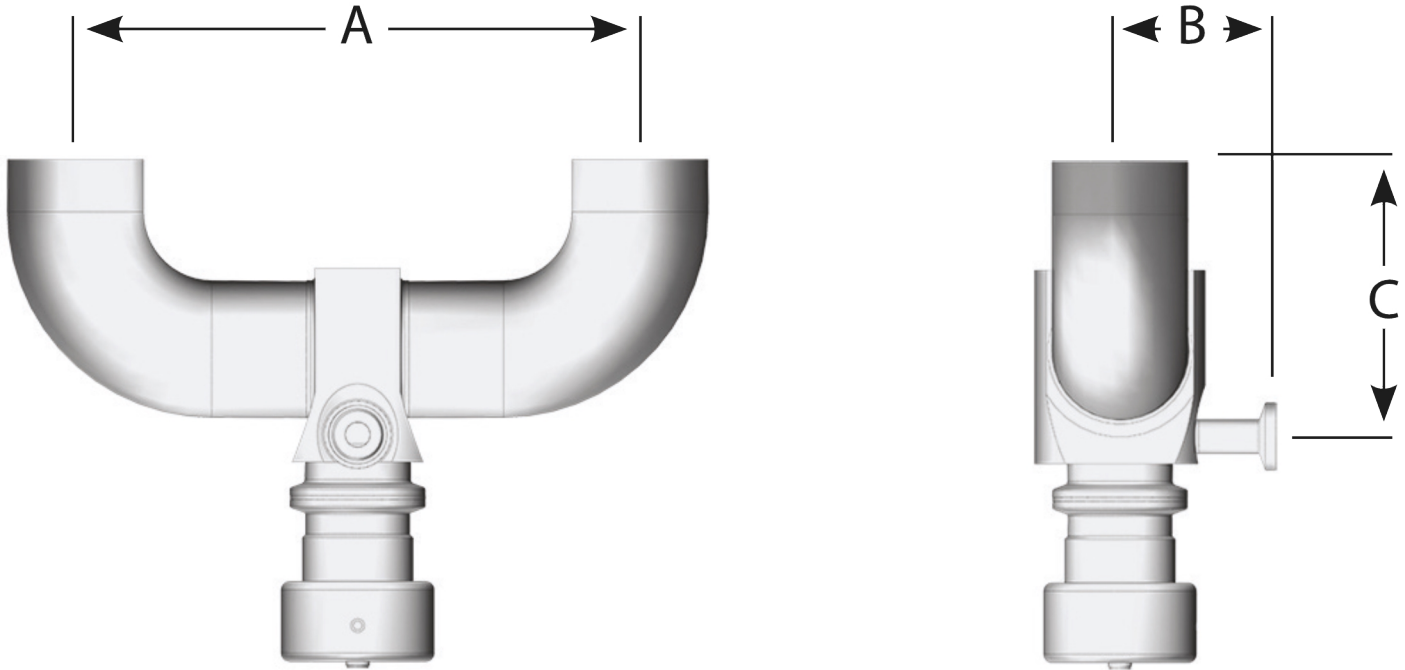

Center to Center Dimension (A)																
Line Size (BW)	0.50		0.75		1.00		1.50		2.00		2.50		3.00		4.00	
Outlet Size (Clamp)	mm	in	mm	in	mm	in	mm	in	mm	in	mm	in	mm	in	mm	in
0.50	195.8	7.71	195.8	7.71	195.8	7.71	223.5	8.80	279.4	11.00	355.6	14.00	393.7	15.50	482.6	19.00
0.75			215.9	8.50	215.9	8.50	254.0	10.00	304.8	12.00	355.6	14.00	393.7	15.50	482.6	19.00
1.00					215.9	8.50	254.0	10.00	304.8	12.00	355.6	14.00	393.7	15.50	482.6	19.00
1.50							303.3	11.94	317.5	12.50	355.6	14.00	406.4	16.00	508.0	20.00
2.00									349.3	13.75	387.4	15.25	406.4	16.00	508.0	20.00
Center to Face Dimension (B)																
Line Size (BW)	0.50		0.75		1.00		1.50		2.00		2.50		3.00		4.00	
Outlet Size (Clamp)	mm	in	mm	in	mm	in	mm	in	mm	in	mm	in	mm	in	mm	in
0.50	52.8	2.08	52.8	2.08	52.8	2.08	52.8	2.08	57.9	2.28	64.3	2.53	64.3	2.53	64.3	2.53
0.75			64.3	2.53	64.3	2.53	64.3	2.53	64.3	2.53	70.6	2.78	70.6	2.78	64.3	2.53
1.00					64.3	2.53	64.3	2.53	64.3	2.53	70.6	2.78	70.6	2.78	64.3	2.53
1.50							85.1	3.35	85.1	3.35	85.1	3.35	85.1	3.35	82.6	3.25
2.00									82.6	3.25	82.6	3.25	82.6	3.25	82.6	3.25
Center to Face Dimension (C)																
Line Size (BW)	0.50		0.75		1.00		1.50		2.00		2.50		3.00		4.00	
Outlet Size (Clamp)	mm	in	mm	in	mm	in	mm	in	mm	in	mm	in	mm	in	mm	in
0.50	89.7	3.53	96.0	3.78	96.0	3.78	121.2	4.77	152.1	5.99	192.5	7.58	217.9	8.58	273.6	10.77
0.75			103.1	4.06	103.1	4.06	131.8	5.19	163.6	6.44	190.5	7.50	220.7	8.69	270.3	10.64
1.00					103.1	4.06	128.8	5.07	160.5	6.32	187.2	7.37	217.7	8.57	267.2	10.52
1.50							147.6	5.81	179.3	7.06	205.0	8.07	230.1	9.06	287.0	11.30
2.00									173.0	6.81	198.6	7.82	223.8	8.81	280.7	11.05

All Process Line and Outlet sizes are in inches

Weights

Weights																
Line Size (BW)	0.50		0.75		1.00		1.50		2.00		2.50		3.00		4.00	
Outlet Size (Clamp)	Kg	lb	Kg	lb	Kg	lb	Kg	lb	Kg	lb	Kg	lb	Kg	lb	Kg	lb
0.50	1.5	3.33	1.8	3.89	2.1	4.63	2.0	4.53	3.0	6.59	4.1	8.99	6.3	13.75	8.6	18.93
0.75			2.8	6.09	2.1	4.63	2.8	6.13	3.8	8.29	4.1	8.99	6.3	13.75	8.6	18.93
1.00					3.9	8.58	3.8	8.28	4.7	10.44	5.1	11.14	7.3	16.00	9.6	21.08
1.50							8.8	19.43	9.7	21.29	10.0	22.09	10.4	22.95	17.0	37.43
2.00									9.5	20.99	9.9	21.79	10.4	22.85	17.0	37.43

All weights are with manual actuator.

All Process Line and Outlet sizes are in inches.

Disclaimer: The information contained in this document is believed to be correct but ASEPCO accepts no liability for any errors it contains and reserves the right to alter specifications without notice. It is the users responsibility to ensure product suitability for use within their application. Radial diaphragm is a trademark of ASEPCO Corporation. Tri-Clamp is a registered trademark of Alfa Laval Corporate AB. A member of Watson-Marlow Fluid Technology Group, A Spirax-Sarco Engineering plc company.



wmfts.com/global
14 July 2025