# 323 Drive with 501RL2 pumphead



300 series cased pump

#### **Features and benefits**

- Manual or remote 4-20 mA or 0-10 V analogue auto-control or RS232 control commands
- Up to 133:1 enhanced digital manual speed control
- Keypad lock prevents tampering or accidental changes
- MemoDose facility for accurate single-shot dosing
- Digital TTL inputs for run/stop and direction
- Instantly reversible
- Two-year warranty



# CE

## **Technical specifications**

	323 Drive with 501RL2 pumphead
Pumphead number of channels	1
Speed control ratio	147:1
Max. operating speed	200 rpm
Min. operating speed	1.5 rpm
Operating temperature range	5 to 40 °C
Operating temperature range	40 to 104 °F
Weight	5.6 kg
Weight	12.3 lbs
Drive control options	U
Drive standards	CE, cETLus
Drive ingress rating	IP31
Drive humidity	80 % up to 31 °C, 88 °F, decreasing linearly to 50 % at 40 °C, 104 °F
Drive noise	<70dB(A) at 1m
Max. altitude	2000 m
Drive power supply	1ph, 100-120 V, 220-240 V 50/60 Hz, 100 VA
Compatible tubing bore size	1.6, 3.2, 4.8, 6.4, 8 mm
Compatible tubing wall thickness	2.4 mm


#### **Materials of construction**

	323 Drive with 501RL2 pumphead			
Drive casework	ABS (Acrylonitrile butadiene styrene), Powder coated aluminium LM24			
Guide roller assembly	Delrin			
Pumphead guard	Polycarbonate			
Pumphead rotor assembly	Mazak, chromium plated			
Spindles	Stainless steel			
Springs	Stainless steel			
Tube clamp assembly	Acetal copolymer (Kematal)			

Information listed covers the complete range.

For detailed specifications of individual models/components refer to user manual or contact WMFTS representative.

# 323 Drive with 501RL2 pumphead dimensions

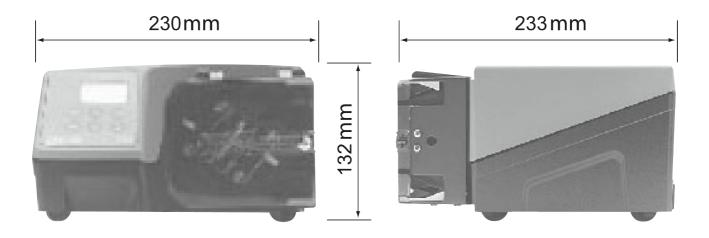

# **Control options**

Functionality	323E	323S	323U	323Du	
Manual control					
Keypad speed control	✓	✓	<b>✓</b>	✓	
15-400rpm 27:1 speed control	✓				
3-400rpm 133:1 speed control		✓	✓	✓	
1.5-220rpm 147:1 speed control		✓	<b>✓</b>	✓	
MemoDose for accurate single shot dosing		✓	<b>✓</b>	✓	
Remote control					
Run/Stop, direction change via contact closure or 5V TTL			<b>✓</b>	✓	
Remote operation of MemoDose (hand/foot switch)			<b>✓</b>	✓	
Analogue speed control					
Inputs; 0-10V or 4-20mA			<b>✓</b>	✓	
RS232 digital communication					
Digital serial control				✓	
Security					
Keypad lock to protect set-up		<b>✓</b>	<b>/</b>	<b>✓</b>	

## **Product codes**

Pumpheads product codes				
Description		Partcode		
501RL2	pumphead for 2.4 mm wall tube	053.0001.L20		

<sup>\*</sup> Replace last 0 with A, E or U for American, European or UK mains supply (American 110 V 1ph 60 Hz)
If required for use with STA-PURE and Chem-Sure tubing, please contact Watson-Marlow Bredel technical support for drive and pumphead selection

Tubing product codes						
Bore/wall (mm)	Bioprene	Marprene	Platinum Silicone			
1.6 / 2.4	933.0016.024	902.0016.024	913.A016.024			
3.2 / 2.4	933.0032.024	902.0032.024	913.A032.024			
4.8 / 2.4	933.0048.024	902.0048.024	913.A048.024			
6.4 / 2.4	933.0064.024	902.0064.024	913.A064.024			
8.0 / 2.4	933.0080.024	902.0080.024	913.A080.024			
Bore/wall (mm)	STA-PURE®	STA-PURE®Series PCS	STA-PURE®Series PFL	Chem-Sure®		
1.6 / 2.4	960.0016.024	961.0016.024	966.0016.024	965.0016.024		
3.2 / 2.4	960.0032.024	961.0032.024	966.0032.024	965.0032.024		
4.8 / 2.4	960.0048.024	961.0048.024	966.0048.024	965.0048.024		
6.4 / 2.4	960.0064.024	961.0064.024	966.0064.024	965.0064.024		
8.0 / 2.4	960.0080.024			965.0080.024		

Disclaimer: All flow rates shown were obtained pumping water at 20 °C (68 °F) with zero suction and delivery heads. The information contained in this document is believed to be correct but Watson-Marlow Limited accepts no liability for any errors it contains and reserves the right to alter specifications without notice. It is the user's responsibility to ensure product suitability for use within their application. Watson-Marlow, LoadSure, Pumpsil, PureWeld XL, Bioprene, Marprene are registered trademarks of Watson-Marlow Limited. Tri-Clamp is a registered trademark of Alfa Laval Corporate AB. GORE and STA-PURE are registered trademarks of W.L. Gore and Associates. Please state the product code when ordering pumps and tubing.

wmfts.com/global



07 August 2023